

NEW ASKIR 230/12V BR

MAIN APPLICATIONS

Ambulance Car

Emergency Dept.


Home-Care

Patient Transfer

Room / Ward

Tracheostomy

Veterinary

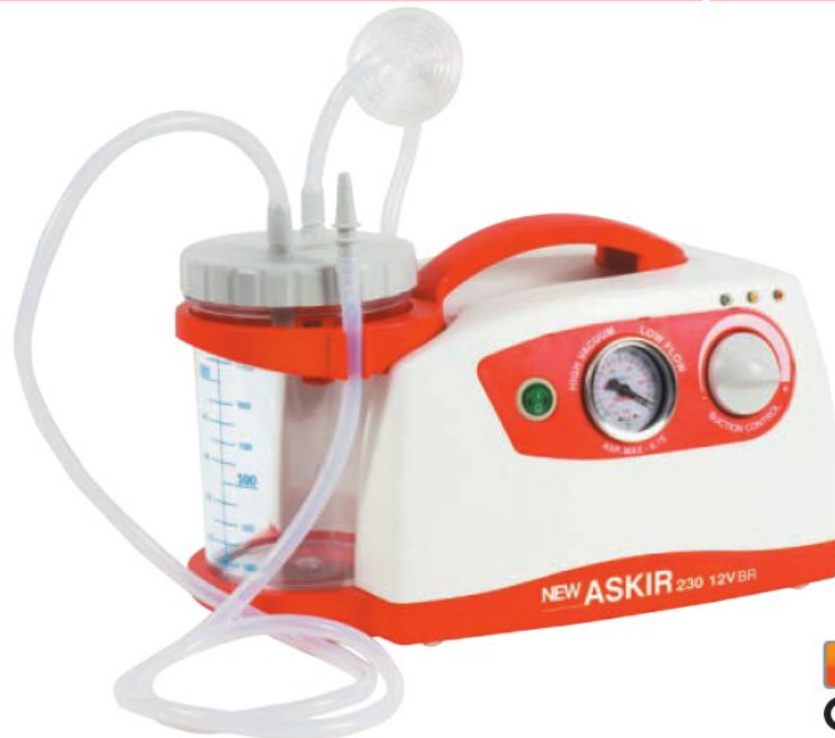

Internal Rechargeable Battery

CE 0123

NEW ASKIR 230/12V BR is an electric suction unit for the aspiration of body liquids, oral, nasal and tracheal aspiration in adults or children. Suitable for emergency applications and for post-operative therapy in home-care. Three different options for operation: AC/DC adapter, rechargeable battery and 12V makes it a versatile suction unit approved for NON-STOP operation without overheating. Provided with visual and acoustic alarm indicating low battery level, the main unit is also equipped with vacuum-meter (bar & kPa) and vacuum regulator. Complies with EN 60601-1-11 "for use in the home healthcare environment" and carries IP21 degree of protection.

ACCESSORIES INCLUDED

- Liquid Collection Jar with overflow valve system (different options, see below)
- 6x10 Silicone Tubes (autoclavable) + Conical Connector \varnothing 8-9-10mm
- Antibacterial & Hydrophobic Filter (single-patient)
- CH20 Canula
- Power Cord with Schuko plug + 12V Car Adapter

OPTIONS

- Stand on Castors REF. 27731
 - Carrying bag* REF. SP 0207
- *for Unit and Consumables

TECHNICAL FEATURES

Motor	Oiless and maintenance-free piston pump		
Power Feeding	230V-50Hz		
	— Internal Rechargeable Battery (12V \approx 2,8 A) - Autonomy 45min - Recharging Time 120÷150min		
	— 12V \approx from Car Battery cigarette lighter		
ISO 10079-1 Classification	HIGH VACUUM / LOW FLOW		
Max Vacuum (adjustable)	-0.75 Bar	-75 kPa	-563 mmHg
Max free air flow rate	16 l/min		
Noise Level	63,0 dB		
Duty cycle	Non-stop operation		
Weight	Kg 4.90		
Size	cm 35 x 21 x 18		
Years of Warranty	2		
Shipping carton	4		
Place of Manufacturing	Italy		

Italian
Medical
Touch
MADE IN ITALY

AVAILABLE CONFIGURATIONS

- RE 310200/01 — NEW ASKIR 230/12V BR with 1000ml autoclavable jar in polycarbonate
- RE 310200/04 — NEW ASKIR 230/12V BR with 2000ml autoclavable jar in polycarbonate
- RE 310201/01 — NEW ASKIR 230/12V BR with 1000ml Flovac[®] disposable liner and reusable container
- RE 310200/25 — NEW ASKIR 230/12V BR with 2000ml Flovac[®] disposable liner and reusable container
 - Jars in APEC[®] autoclavable up to 143°C available upon request